

NBP1-58872**SH3BP2 Antibody****Unit Size**

0.05 mg

Concentration

Lyoph

Storage

Store at -20C. Avoid freeze-thaw cycles.

Applications

WB, IHC-P

Species

Human, Mouse, Rat

Host

Rabbit

Clonality

Polyclonal

Purity

Peptide affinity purified

Preservative

No Preservative

Immunogen

Synthetic peptides corresponding to SH3BP2(SH3-domain binding protein 2) The peptide sequence was selected from the middle region of SH3BP2. Peptide sequence RQPSQADTGGDDSEDEYKVLPLNSVVFVNTTESCEVERLFKATSPRGEPQ.

Recommended Dilutions

Immunohistochemistry-Paraffin, Western Blot

Buffer

PBS & 2% Sucrose lyophilized with the antibody. Centrifuge the vial of lyoph antibody at 12,000 x g for 20 seconds. Add 50ul of distilled water. Vortex followed by centrifuge again to pellet the solution. Final concentration is 1 mg/ml in PBS buffer

Uses


This is a rabbit polyclonal antibody against SH3BP2 and was validated on Western blot.

Images

Western Blot: SH3BP2 Antibody [NBP1-58872] -
 Titration: 1 ug/ml Positive Control: 721_B cell lysate.



Immunohistochemistry-Paraffin: SH3BP2 Antibody [NBP1-58872] - Human kidney Tissue, antibody concentration 5 ug/ml.

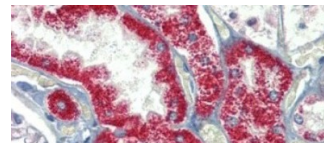

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

