

NBP1-50777

SPHK2 Antibody

Unit Size

0.1 mg

Concentration

0.5 mg/ml

Storage

Aliquot and store at -20C or -80C. Avoid freeze-thaw cycles.

Applications

WB

Species

Human, Mouse, Rat

Host

Rabbit

Clonality

Polyclonal

Purity

Affinity purified

Preservative

0.01% Thimerosal

Immunogen

Synthetic peptide corresponding to residues surrounding amino acid 603 of rat Sphingosine Kinase 2

Recommended Dilutions

Western Blot

Buffer


phosphate-buffered saline (PBS) containing 30% glycerol, 0.5% BSA and 0.01% Thimerosal.

Uses

Western Blot (0.5-4 ug/ml). The antibody recognizes ~69 kDa bands from samples of human, mouse and rat origins.

Images

Western Blot: SPHK2 Antibody [NBP1-50777] - Jurkat cell lysate.

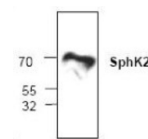

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

