

NBP1-41202**RbAp48 Antibody (15G12)****Unit Size**

0.1 mg

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Store at 4°C short term. Aliquot and store at -20°C long term. Avoid freeze-thaw cycles.

Applications

WB, ICC/IF, IHC, IP

Species

Human

Host

Mouse

Clonality

Monoclonal

Isotype

IgG

Purity

Protein G purified

Preservative

No Preservative

Specificity/Sensitivity

This antibody will recognize human RbAp48 + RbAp46.

Immunogen

The complete coding region (amino acids 1-425) of RbAp48 expressed in E. coli.

Recommended Dilutions

Immunocytochemistry/Immunofluorescence, Immunohistochemistry, Immunoprecipitation, Western Blot

Buffer


Phosphate-buffered saline, pH 7.4, containing no preservatives

Uses

For IF: Use at a concentration of 0.2 - 1 µg/ml. For IP: Use at a concentration of 0.2 - 1 µg/ml. For WB: Use at a concentration of 0.2 - 1 µg/ml. Predicted molecular weight: 47.6 kDa. Not tested in other applications. Optimal dilutions/concentrations should be determined by the researcher.

Images (more available at www.novusbio.com/NBP1-41202)

Western Blot: RbAp48 Antibody (15G12) [NBP1-41202] - Sample (30 µg of whole cell lysate) A: Raji B: K562 C: THP-1 D: HL-60 10% SDS PAGE; antibody diluted at 1:1000.

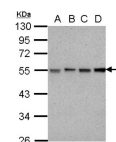

Immunohistochemistry: RbAp48 Antibody (15G12) [NBP1-41202] - Analysis of paraffin-embedded C2C12 xenograft, using RbAp48 antibody at 1:200 dilution.



Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

