

NBP1-32499**SLC25A33 Antibody****Unit Size**

0.1 mg

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Aliquot and store at -20C or -80C. Avoid freeze-thaw cycles.

Applications

WB, IHC, IHC-P

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.01% Thimerosal

Immunogen

Recombinant fragment corresponding to a region within amino acids 114 and 321 of SLC25A33 (Uniprot ID#Q9BSK2)

Recommended Dilutions

Immunohistochemistry, Immunohistochemistry-Paraffin 1:100-1:1000, Western Blot 1:500-1:3000

Buffer

0.1M Tris, 0.1M Glycine, 10% Glycerol (pH7.0).

Uses


For IHC-Paraffin, suggested antigen retrieval using heat mediated 10mM Citrate buffer (pH6.0) or Tris-EDTA buffer (pH8.0).

References


Franzolin E, Miazzi C, Frangini M et al. The pyrimidine nucleotide carrier PNC1 and mitochondrial trafficking of thymidine phosphates in cultured human cells. Exp Cell Res 2012 Oct [PMID: 22677043]

Images

Western Blot: SLC25A33 Antibody [NBP1-32499] - Sample(30 ug whole cell lysate)A: 293TB: A431 G: Raji 10% SDS PAGE, antibody diluted at 1:1000.

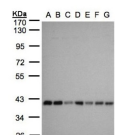

Immunohistochemistry-Paraffin: SLC25A33 Antibody [NBP1-32499] - Paraffin-embedded SAOS-2 xenograft, using antibody at 1:500 dilution.

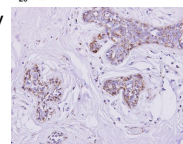

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

