

NBP1-32676**UGP2 Antibody****Unit Size**

0.1 ml

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Store at -20C. Avoid freeze-thaw cycles.

Applications

WB, ICC/IF, IHC, IHC-P

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Affinity purified

Preservative

0.01% Thimerosal

Immunogen

Recombinant fragment corresponding to a region within amino acids 219 and 508 of UGP2 (Uniprot ID#Q16851)

Recommended Dilutions

Immunocytochemistry/Immunofluorescence 1:100-1:200, Immunohistochemistry, Immunohistochemistry-Paraffin 1:100-1:500, Western Blot 1:500-1:3000

Buffer


0.1M Tris, 0.1M Glycine, 10% Glycerol (pH7.0).

Uses

Western blot. The usefulness of this product in other applications has not been determined.

Images (more available at www.novusbio.com/NBP1-32676)

Western Blot: UGP2 Antibody [NBP1-32676] -
 Sample (30 ug of whole cell lysate)A: HeLa S37C 5%
 SDS PAGE, antibody diluted at 1:1000.

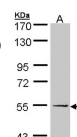

Immunocytochemistry/Immunofluorescence: UGP2
 Antibody [NBP1-32676] - Paraformaldehyde-fixed
 HeLa, using antibody at 1:200 dilution.



Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

