

NBP1-32791**O-GlcNAc Transferase p110 subunit Antibody****Unit Size**

0.1 mg

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Store at -20C. Avoid freeze-thaw cycles.

Applications

WB, ICC/IF, IHC-P

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Affinity purified

Preservative

0.01% Thimerosal

Marker

Golgi Apparatus Marker

Immunogen

Recombinant fragment corresponding to a region within amino acids 223 and 502 of O-GlcNAc transferase (Uniprot ID#O15294)

Recommended Dilutions

Immunocytochemistry/Immunofluorescence, Immunohistochemistry-Paraffin, Western Blot 1:500-1:3000

Buffer

0.1M Tris, 0.1M Glycine, 10% Glycerol (pH7.0).

Uses


Western blot. The usefulness of this product in other applications has not been determined.

References

Ise H, Goto M, Komura K et al. Engulfment and clearance of apoptotic cells based on a GlcNAc-binding lectin-like property of surface vimentin. Glycobiology 2012 Jun [PMID: 22345628] (WB)

Images (more available at www.novusbio.com/NBP1-32791)

Western Blot: O-GlcNAc Transferase p110 subunit Antibody [NBP1-32791] - Sample (30 ug of whole cell lysate)A: HeLa 7.5% SDS PAGE, antibody diluted at 1:1000.

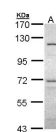

Immunocytochemistry/Immunofluorescence: O-GlcNAc Transferase p110 subunit Antibody [NBP1-32791] - Paraformaldehyde-fixed HeLa, using antibody at 1:200 dilution.



Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA
 8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com
Novus Canada
 461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com
Novus Europe
 12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com
