

NBP1-31540**Nuclear Factor Erythroid Derived 2 Antibody****Unit Size**

0.1 ml

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Store at -20C. Avoid freeze-thaw cycles.

Applications

WB

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Affinity purified

Preservative

0.01% Thimerosal

Immunogen

Recombinant fragment corresponding to a region within amino acids 1 and 373 of NF-E2 (Uniprot ID#Q16621)

Recommended Dilutions

Western Blot 1:500-1:3000

Buffer


0.1M Tris, 0.1M Glycine, 10% Glycerol (pH7.0).

Uses


Western blot

Images

Western Blot: Nuclear Factor Erythroid Derived 2 Antibody [NBP1-31540] - Sample (30 ug of whole cell lysate) A: 293T 10% SDS PAGE; antibody diluted at 1:1000.

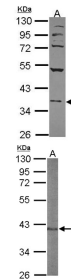

Western Blot: Nuclear Factor Erythroid Derived 2 Antibody [NBP1-31540] - Sample (50 ug of whole cell lysate) A: mouse spleen 10% SDS PAGE; antibody diluted at 1:1000.

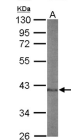

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

