

NB100-91878**MMP3 Antibody****Unit Size**

0.1 mg

Concentration

1.0 mg/ml

Storage

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

Applications

WB, ICC/IF, IHC, IHC-P

Species

Human, Mouse, Rat

Molecular Weight

53kDa

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.05% Sodium Azide

Immunogen

Synthetic peptide corresponding to the residues surrounding Serine 453 of human MMP3.

Recommended Dilutions

Immunocytochemistry/Immunofluorescence 1:50-1:200, Immunohistochemistry, Immunohistochemistry-Paraffin 1:50-1:200, Western Blot 1:500- 1:1000

Buffer


Phosphate buffered saline, pH 7.2.

Uses


Western blot (WB) analysis of MMP-3 (S453) pAb in extracts from 293 cells. Immunohistochemistry (IHC) analyzes of MMP-3 (S453) pAb in paraffin-embedded human breast carcinoma tissue.

Images

Western Blot: MMP3 Antibody [NB100-91878] - Analysis of MMP-3 (S453) antibody in extracts from 293 cells.

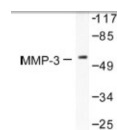

Immunohistochemistry: MMP3 Antibody [NB100-91878] - Analyzes of MMP-3 (S453) antibody in paraffin-embedded human breast carcinoma tissue.

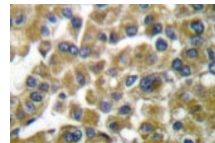

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

