

NB100-55328**SMCX Antibody****Unit Size**

0.1 ml

Concentration

0.2 mg/ml

Storage

Store at 4°C. Do not freeze.

Applications

WB, IP

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.09% Sodium Azide

Immunogen

The immunogen recognized by this antibody maps to a region between residue 1510 and the C-terminus (residue 1560) of human Jumonji, AT rich interactive domain 1C using the numbering given in entry NP_004178.2 (GeneID 8242).

Recommended Dilutions

Immunoprecipitation 2 to 5 mcg/mg lysate, Western Blot 1:2000 to 1:10000

Buffer

Tris-buffered Saline and 0.1% BSA

Uses


This antibody is useful for Western Blot and Immunoprecipitation

References

M. L. C. Powell, J. A. Smith, M. E. Sowa, et al. NCoR1 Mediates Papillomavirus E8^{E2C} Transcriptional Repression. J Virol ; 84 (9): 4451-4460. 2010. [PMID: 20181716]

Images

Western Blot: SMCX Antibody [NB100-55328] - Detection of Human JARID1C on HeLa whole cell lysate using NB100-55328. JARID1C was also immunoprecipitated by rabbit anti-JARID1C antibody NB100-55327.

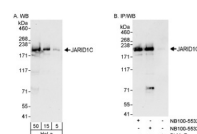

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

