

NB100-59846**RBM12 Antibody****Unit Size**

0.1 ml

Concentration

0.2 mg/ml

Storage

Store at 4°C. Do not freeze.

Applications

WB, IHC, IHC-P, IP

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.09% Sodium Azide

Immunogen

The immunogen recognized by this antibody maps to a region between residue 625 and 675 of human RNA binding motif protein 12 using the numbering given in entry NP_690051.1 (GeneID 10137).

Recommended Dilutions

Immunohistochemistry, Immunohistochemistry-Paraffin 1:100 - 1:500, Immunoprecipitation 2 to 5 ug/mg lysate, Western Blot 1:2000 to 1:10000

Buffer


Tris-buffered Saline and 0.1% BSA

Uses

Western Blot and Immunoprecipitation

Images (more available at www.novusbio.com/NB100-59846)

Western Blot: RBM12 Antibody [NB100-59846] - Detection of Human RBM12 on HeLa whole cell lysate using NB100-59846. RBM12 was also immunoprecipitated by rabbit anti-RBM12 antibody NB100-59847.

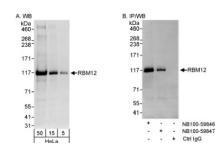

Immunohistochemistry-Paraffin: RBM12 Antibody [NB100-59846] - FFPE section of human prostate carcinoma. Affinity purified rabbit anti-RBM12 used at a dilution of 1:250.



Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

