

NB100-1565**AF9 Antibody****Unit Size**

0.1 mg

Concentration

1.0 mg/ml

Storage

Store at 4°C. Do not freeze.

Applications

WB, IHC-P, IP

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.09% Sodium Azide

Immunogen

The immunogen recognized by this antibody maps to a region between residues 250 and 300 of human ALL1 Fused Gene from Chromosome 9 Protein using the numbering given in entry NP_004520.1 (GeneID 4300).

Recommended Dilutions

Immunoprecipitation 1 -4 ug/mg of lysate, Western Blot 1:5000 - 1:15000, Immunohistochemistry-Paraffin

Buffer

Tris-citrate/phosphate buffer, pH 7 to 8


Uses

This antibody is useful for Western blot and Immunoprecipitation

Images

Western Blot: AF9 Antibody [NB100-1565] - Detection of Human AF9 on HeLa whole cell lysates using NB100-1565. AF9 was also immunoprecipitated using rabbit anti-AF9 antibodies NB100-1564, NB100-1566 and another AF9 Ab.

Immunohistochemistry-Paraffin: AF9 Antibody [NB100-1565] - Human ovarian carcinoma. Antibody used at a dilution of 1:1000 (1ug/ml).

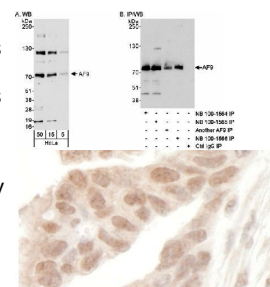

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

