

NB100-314**AIB1 Antibody****Unit Size**

0.1 ml

Concentration

Please see the vial label for concentration.

Storage

Store at 4°C. Do not freeze.

Applications

WB, IP, PLA

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.1% Sodium Azide

Immunogen

A synthetic peptide, which represented a portion of human nuclear receptor coactivator 3 encoded within exon 13 (LocusLink ID 8202).

Recommended Dilutions

Immunoprecipitation 1 -4 ug/mg lysate, Proximity Ligation Assay 1:1000 -1:20000, Western Blot 1:1000 -1:10000

Buffer


Tris-citrate/phosphate buffer, pH 7 to 8

Uses

NB 100-314 may be used for Western blot and immunoprecipitation. This protein is a cytoplasmic coactivator that shifts to the nucleus upon activation. It can be immunoprecipitated and analyzed on Western blot using nuclear or cytosolic extracts. Suggested working dilutions: * Western Blot - 1:1,000 ~ 1:10,000 (5 ug nuclear extract) Immunocytochemistry - Not Done Immunoprecipitation - 1-4 ug/mg ly... See more at www.novusbio.com/NB100-314

Images

Western Blot: AIB1 Antibody [NB100-314] - Nuclear extract (5 mcg for WB, 10 mg for CoIP) from HeLa cells. Affinity purified rabbit anti-SRC3 used at 0.33 mcg/ml for WB and 20 mcg/10 mg extract for CoIP.

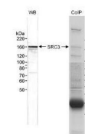

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

