

Recombinant hSUMO2-AMC


CATALOG #:	6413-50	50 µg
ALTERNATE NAMES:	Small ubiquitin-related modifier 2, SUMO-2, Ubiquitin-like protein SMT3B, SMT3 homolog 2, Sentrin-2, HSMT3, SUMO-3, SUMO2, SMT3B, SMT3H2, MGC117191.	
PURITY:	≥ 95% by RP-HPLC	
MOL. WEIGHT:	10.619 kDa	
FORMULATION:	Lyophilized powder.	
STORAGE CONDITIONS:	Aliquot and store at -80°C. Avoid repeated freezing and thawing cycles.	

DESCRIPTION:

hSUMO2-AMC is prepared via the conjugation of 7-amino-4-methylcoumarin (AMC) to the C-terminal Gly of hSUMO2. This conjugation quenches the intrinsic fluorescence of AMC. Upon incubation with a protease recognizing hSUMOs, such as Senp1 or Senp2, AMC is released and the increase in fluorescence at 460 nm (Exc. 380 nm) can be measured. The enhanced purity of hSUMO2-AMC compared to that of other commercial vendors leads to a greater signal to noise ratio. This protein contains no extraneous tags.

APPLICATIONS:

- Used for determining the activity and specificity of SUMO-specific hydrolases. High throughput screening for modulators of desumoylase activity.

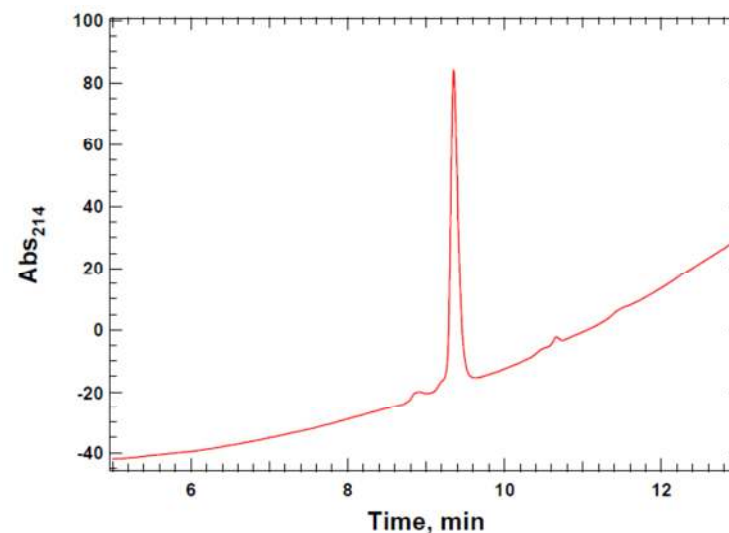

RP-HPLC of Human Recombinant hSUMO2-AMC

RELATED PRODUCTS:

- SUMO1, human recombinant (Cat. No. 4941-100, -1000)
- SUMO2, human recombinant (Cat. No. 4942-100, -1000)
- SUMO3, human recombinant (Cat. No. 4943-100, -1000)
- hSUMO1-AMC (Cat. No. 6412-50)
- hSUMO2-Rhodamine (Cat. No. 6414-50)

FOR RESEARCH USE ONLY! Not to be used in humans.