# *BioVision* ACAA2 Polyclonal Antibody

CATALOG #:	3538-100
AMOUNT:	100 µg
LOT #:	
IMMUNOGEN:	Synthetic peptide surrounding amino acid 205 of rat ACAA2
INTERNAL ID#:	BV-K67
HOST:	Rabbit
SPECIES REACTIVITY:	human, mouse, rat

## FORMULATION:

100 µg (0.5 mg/ml) affinity purified rabbit anti-ACAA2 polyclonal antibody in phosphate buffered saline (PBS), pH 7.2, containing 30% glycerol, 0.5% BSA, 5 mM EDTA and 0.01% thimerosal.

## STORAGE CONDITIONS:


Store at -20°C. For long term storage, aliquot and freeze at -70°C. Avoid repeated freeze/defrost cycles.

### DESCRIPTION:

ACAA2 (acetylcoenzyme A acyltransferase 2), is member of the thiolase family of enzymes and is involved in lipid metabolism. ACAA2 is localized in the mitochondria. ACCA2 catalyzes the last conversion step of acetyl-CoA to 3-oxoacyl-CoA, in the fatty acid oxidation pathway.

## APPLICATION AND USAGE:

The antibody can be used in Western Blot analysis (1-4 µg/ml). However, the optimal concentrations should be determined individually. Blocking peptide (Cat.# 3538BP-50) is available separately.

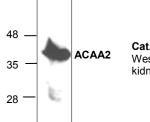

Cat. #: 3538-100 ACAA2 Antibody. Western blot analysis of ACAA2 using rat kidney tissue lysate.

FOR RESEARCH USE ONLY! Not to be used in humans.

### **RELATED PRODUCTS:**

Apoptosis Detection Kits & Reagents

- Annexin V Kits & Bulk Reagents
- Caspase Assay Kits & Reagents
- Mitochondrial Apoptosis Kits & Reagents
- Nuclear Apoptosis Kits & Reagents
- Apoptosis Inducers and Set
- Apoptosis siRNA Vectors

Cell Fractionation System

- Mitochondria/Cytosol Fractionation Kit
- Nuclear/Cytosol Fractionation Kit
- Membrane Protein Extraction Kit
- Cytosol/Particulate Rapid Separation Kit
- Mammalian Cell Extraction Kit
- FractionPREP Fractionation System
- Cell Proliferation & Senescence
  - Quick Cell Proliferation Assay Kit
  - Senescence Detection Kit
  - High Throughput Apoptosis/Cell Viability Assay Kits
  - LDH-Cytotoxicity Assay Kit
  - Bioluminescence Cytotoxicity Assay Kit
  - Live/Dead Cell Staining Kit

Cell Damage & Repair

- HDAC Fluorometric & Colorimetric Assays & Drug Discovery Kits
- HAT Colorimetric Assay Kit & Reagents
- DNA Damage Quantification Kit
- Glutathione & Nitric Oxide Fluorometric & Colorimetric Assay Kits

Signal Transduction

- cAMP & cGMP Assay Kits
- Akt & JNK Activity Assay Kits
- Beta-Secretase Activity Assay Kit

Adipocyte & Lipid Transfer

- Recombinant Adiponectin, Survivin, & Leptin
- CETP Activity Assay & Drug Discovery Kits
- PLTP Activity Assay & Drug Discovery Kits
- Total Cholesterol Quantification Kit
- Molecular Biology & Reporter Assays
  - siRNA Vectors
  - Cloning Insert Quick Screening Kit
  - Mitochondrial & Genomic DNA Isolation Kits
  - 5 Minutes DNA Ligation Kit
  - 20 Minutes Gel Staining/Destaining Kit
  - β-Galactosidase Staining Kit & Luciferase Reporter Assay Kit

Growth Factors and Cytokines Monoclonal and Polyclonal Antibodies