

Bid Polyclonal Antibody

CATALOG #: 3272-100

AMOUNT: 100 µg

LOT #: _____

HOST: Rabbit

IMMUNOGEN: Synthetic peptide surrounding amino acid 45 of mouse Bid (ID# BV-550).

SPECIES REACTIVITY: Mouse

FORMULATION:

100 µg (0.5 mg/ml) affinity purified rabbit anti-Bid polyclonal antibody in phosphate buffered saline (PBS), pH 7.2, containing 30% glycerol, 0.5% BSA, 0.01% thimerosal.

STORAGE CONDITIONS:

Store at -20°C. For long term storage, aliquot and freeze at -70°C. Avoid repeated freeze/thaw cycles.

BACKGROUND DESCRIPTION:


Bid, a BH3 domain-containing proapoptotic Bcl-2 family member, is localized in the cytosolic fraction of cells as an inactive precursor. Its active form is generated upon proteolytic cleavage by caspase-8 in the Fas signaling pathway. Cleaved Bid translocates to mitochondria and releases its potent proapoptotic activity, which in turn induces cytochrome c release and mitochondrial damage. The cytochrome c releasing activity of Bid was antagonized by Bcl-2. Mutation in the SH3 domain can diminish the cytochrome c releasing activity. In the animal model studies, Bid-deficient mice are found resistant to the lethal effects of death factor signals relayed through Fas.

SPECIFICITY:

The antibody detects the full length (22 kDa) mouse Bid. Reactivity to other species has not been tested.

APPLICATION AND USAGE:

The antibody can be used in Western blot analysis (0.5-5 µg/ml). Per researcher's feedback, it also works well in immunoprecipitation (20 µg/ml) and immunohistochemistry (20 µg/ml). However, the optimal conditions should be determined individually. Blocking peptide (Cat.# 3272BP-50) is available separately.

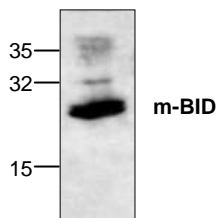

Cat. #: 3272-100
Western blot analysis of BID expression in 3T3 cell lysate.

RELATED PRODUCTS:

Apoptosis Detection Kits & Reagents

- Annexin V Kits & Bulk Reagents
- Caspase Assay Kits & Reagents
- Mitochondrial Apoptosis Kits & Reagents
- Nuclear Apoptosis Kits & Reagents
- Apoptosis Inducers and Set
- Apoptosis siRNA Vectors

Cell Fractionation System

- Mitochondria/Cytosol Fractionation Kit
- Nuclear/Cytosol Fractionation Kit
- Membrane Protein Extraction Kit
- Cytosol/Particulate Rapid Separation Kit
- Mammalian Cell Extraction Kit
- FractionPREP Fractionation System

Cell Proliferation & Senescence

- Quick Cell Proliferation Assay Kit
- Senescence Detection Kit
- High Throughput Apoptosis/Cell Viability Assay Kits
- LDH-Cytotoxicity Assay Kit
- Bioluminescence Cytotoxicity Assay Kit
- Live/Dead Cell Staining Kit

Cell Damage & Repair

- HDAC Fluorometric & Colorimetric Assays & Drug Discovery Kits
- HAT Colorimetric Assay Kit & Reagents
- DNA Damage Quantification Kit
- Glutathione & Nitric Oxide Fluorometric & Colorimetric Assay Kits

Signal Transduction

- cAMP & cGMP Assay Kits
- Akt & JNK Activity Assay Kits
- Beta-Secretase Activity Assay Kit

Adipocyte & Lipid Transfer

- Recombinant Adiponectin, Survivin, & Leptin
- CETP Activity Assay & Drug Discovery Kits
- PLTP Activity Assay & Drug Discovery Kits
- Total Cholesterol Quantification Kit

Molecular Biology & Reporter Assays

- siRNA Vectors
- Cloning Insert Quick Screening Kit
- Mitochondrial & Genomic DNA Isolation Kits
- 5 Minutes DNA Ligation Kit
- 20 Minutes Gel Staining/Destaining Kit
- β-Galactosidase Staining Kit & Luciferase Reporter Assay Kit

Growth Factors and Cytokines

Monoclonal and Polyclonal Antibodies

FOR RESEARCH USE ONLY! Not to be used on humans.