

SGK Polyclonal Antibody

CATALOG #:	3459-100
AMOUNT:	100 µg
LOT #:	_____
HOST:	Rabbit
IMMUNOGEN:	Synthetic peptide corresponding to the C-terminus of human SGK (Internal ID# BV-579).
SPECIES REACTIVITY:	Human, Mouse, Rat

FORMULATION:

100 µg (0.5 mg/ml) affinity purified rabbit anti-SGK polyclonal antibody in phosphate buffered saline (PBS), pH 7.2, containing 30% glycerol, 0.5% BSA, 0.01% thimerosal.

STORAGE CONDITIONS:

Store at -20°C. For long-term storage, aliquot and freeze at -70°C. Avoid repeated freeze/thaw cycles.

BACKGROUND:


Serum and glucocorticoid-inducible kinase (SGK) is a serine/threonine kinase that is closely related to Akt. SGK is rapidly induced in response to a variety of stimuli, including serum, hormones, and osmotic shock, etc. SGK has been involved in modulating anti-apoptotic and cell cycle regulation, as well as in activating certain potassium, sodium and chloride channels.

SPECIFICITY:

The antibody detects SGK from samples of human, mouse, and rat origins.

APPLICATION AND USAGE:

The antibody can be used for Western blotting (0.5-4 µg/ml). Per researcher's feedback, it can also be used in immunoprecipitation (10-20 µg/ml). However, the optimal conditions should be determined individually. Blocking peptide (Cat.# 3459BP-50) is available separately.

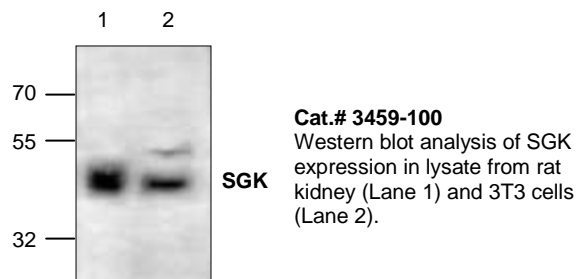

FOR RESEARCH USE ONLY! Not to be used on humans.

RELATED PRODUCTS:

Apoptosis Detection Kits & Reagents

- Annexin V Kits & Bulk Reagents
- Caspase Assay Kits & Reagents
- Mitochondrial Apoptosis Kits & Reagents
- Nuclear Apoptosis Kits & Reagents
- Apoptosis Inducers and Set
- Apoptosis siRNA Vectors

Cell Fractionation System

- Mitochondria/Cytosol Fractionation Kit
- Nuclear/Cytosol Fractionation Kit
- Membrane Protein Extraction Kit
- Cytosol/Particulate Rapid Separation Kit
- Mammalian Cell Extraction Kit
- FractionPREP Fractionation System

Cell Proliferation & Senescence

- Quick Cell Proliferation Assay Kit
- Senescence Detection Kit
- High Throughput Apoptosis/Cell Viability Assay Kits
- LDH-Cytotoxicity Assay Kit
- Bioluminescence Cytotoxicity Assay Kit
- Live/Dead Cell Staining Kit

Cell Damage & Repair

- HDAC Fluorometric & Colorimetric Assays & Drug Discovery Kits
- HAT Colorimetric Assay Kit & Reagents
- DNA Damage Quantification Kit
- Glutathione & Nitric Oxide Fluorometric & Colorimetric Assay Kits

Signal Transduction

- cAMP & cGMP Assay Kits
- Akt & JNK Activity Assay Kits
- Beta-Secretase Activity Assay Kit

Adipocyte & Lipid Transfer

- Recombinant Adiponectin, Survivin, & Leptin
- CETP Activity Assay & Drug Discovery Kits
- PLTP Activity Assay & Drug Discovery Kits
- Total Cholesterol Quantification Kit

Molecular Biology & Reporter Assays

- siRNA Vectors
- Cloning Insert Quick Screening Kit
- Mitochondrial & Genomic DNA Isolation Kits
- 5 Minutes DNA Ligation Kit
- 20 Minutes Gel Staining/Destaining Kit
- β -Galactosidase Staining Kit & Luciferase Reporter Assay Kit

Growth Factors and Cytokines

Monoclonal and Polyclonal Antibodies