


LYZL6 Antibody (Center)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)

Catalog #	Applications:	Reactivity:	Accessions:
AP16919c	WB, E	H	O75951 , NP_065159.1
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB36229

Application Data: Calculated MW: 16956 Da

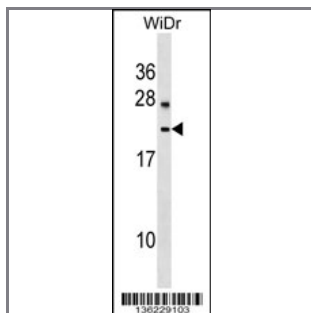

LYZL6 Antibody (Center) (Cat. #AP16919c) western blot analysis in WiDr cell line lysates (35ug/lane). This demonstrates the LYZL6 antibody detected the LYZL6 protein (arrow).

Gene ID:	Gene Symbol:
57151	LYZL6

Other Names:

Lysozyme-like protein 6;LYZL6;LYC1

Target/Specificity:

This LYZL6 antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 36-65 amino acids from the Central region of human LYZL6.

Application Notes:

The suggested dilution is:

ELISA 1:1,000

Western blotting 1:100~500

Format:

Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.

Storage:

LYZL6 Antibody (Center) can be refrigerated at 2-8°C for up to 6 months. For long term storage,

place the at -20°C.

Precautions:

LYZL6 Antibody (Center) is for research use only and not for use in diagnostic or therapeutic procedures.
