

NBP1-71890**MARK2 Antibody****Unit Size**

0.1 mg

Concentration

1 mg/ml

Storage

Store at 4°C. Do not freeze.

Applications

WB, IP

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.09% Sodium Azide

Immunogen

The immunogen for this product maps to a region between residue 350 and 400 of human MAP/Microtubule Affinity-Regulating Kinase 2 using the numbering given in entry NP_001034558.2 (GeneID 2011).

Recommended Dilutions

Immunoprecipitation 2 - 10 mcg/mg, Western Blot 1:2000 - 1:10000

Buffer


Tris-citrate/phosphate buffer, pH 7 to 8

Uses

This antibody works in Western Blot, Immunoprecipitation.

Images

Western Blot: MARK2 Antibody [NBP1-71890] - Whole cell lysate from 293T (5, 15 and 50 mcg for WB; 1 mg for IP, 20% of IP loaded) and Jurkat (J; 50 mcg) cells. Affinity purified rabbit anti-MARK2 antibody used for WB at 0.1 mcg/ml (A) and 1 ... See more at www.novusbio.com/NBP1-71890

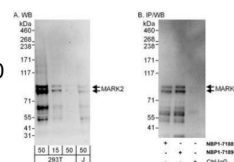

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

