

NBP1-33549**SSB Antibody****Unit Size**

0.1 ml

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Store at -20C. Avoid freeze-thaw cycles.

Applications

WB, ICC/IF, IHC, IHC-P

Species

Human, Mouse

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.01% Thimerosal

Immunogen

Recombinant fragment corresponding to a region within amino acids 6 and 242 of SSB.

Recommended Dilutions

Immunocytochemistry/Immunofluorescence 1:100-1:1000, Immunohistochemistry 1:100-1:1000, Immunohistochemistry-Paraffin 1:100-1:1000, Western Blot 1:1000-1:10000

Buffer

0.1M Tris, 0.1M Glycine, 20% Glycerol (pH7).

Uses


For IHC: Use suggested antigen retrieval using heat mediated 10mM Citrate buffer (pH6.0) or Tris-EDTA buffer (pH8.0).

References

Wang Q, Liu Y, An D et al. Regulation of hepatitis C virus translation initiation by iron: role of eIF3 and La protein. Virus Res 2012 Aug [PMID: 22634302]

Images (more available at www.novusbio.com/NBP1-33549)

Western Blot: SSB Antibody [NBP1-33549] - Sample (30 ug of whole cell lysate) A: 293T 10% SDS PAGE; antibody diluted at 1:1000.

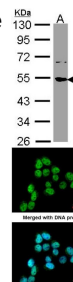

Immunocytochemistry/Immunofluorescence: SSB Antibody [NBP1-33549] - Analysis of paraformaldehyde-fixed A431, using SSB (NBP1-33549) antibody at 1:200 dilution.

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA
 8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com
Novus Canada
 461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com
Novus Europe
 12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com
