

NBP1-31254**DUSP26 Antibody****Unit Size**

0.1 ml

Concentration

Concentrations vary lot to lot. See vial label for concentration.

Storage

Store at -20C. Avoid freeze-thaw cycles.

Applications

WB, DB, ICC/IF

Species

Human, Mouse

Molecular Weight

24kDa

Host

Rabbit

Clonality

Polyclonal

Purity

Affinity purified

Preservative

0.01% Thimerosal

Immunogen

Synthetic peptide corresponding to a region within amino acids 22 and 117 of DUSP26 (Uniprot ID#Q9BV47)

Recommended Dilutions

Dot Blot, Immunocytochemistry/Immunofluorescence, Western Blot 1:500-1:3000

Buffer


0.1M Tris, 0.1M Glycine, 10% Glycerol (pH7.0).

Uses


Dot blot, Western blot. The usefulness of this product in other applications has not been determined.

Images (more available at www.novusbio.com/NBP1-31254)

Western Blot: DUSP26 Antibody [NBP1-31254] -
 Sample (50 ug of whole cell lysate) A: mouse brain
 12% SDS PAGE, antibody diluted at 1:500.

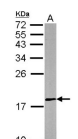

Immunocytochemistry/Immunofluorescence:
 DUSP26 Antibody [NBP1-31254] -
 Paraformaldehyde-fixed HeLa, using antibody at
 1:200 dilution.

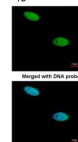

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

