

**NBP1-21372****TAF2 Antibody****Unit Size**

0.1 mg

**Concentration**

1.0 mg/ml

**Storage**

Store at 4C. Do not freeze.

**Applications**

IP, WB (-)

**Species**

Human

**Host**

Rabbit

**Clonality**

Polyclonal

**Purity**

Immunogen affinity purified

**Preservative**

0.09% Sodium Azide

**Immunogen**

The immunogen recognized by this antibody maps to a region between residue 1149 and 1199 of human TBP-associated factor 2 using the numbering given in entry NP\_003175.1 (GeneID 6873).

**Recommended Dilutions**

Immunoprecipitation 5 -10 ug/mg lysate  
 Western Blot (Negative)

**Buffer**


Tris-citrate/phosphate buffer, pH 7 to 8

**Uses**

This antibody is useful for immunoprecipitation. It is not recommended for western analysis.

**Images**

Immunoprecipitation: TAF2 Antibody [NBP1-21372] - Whole cell lysate (1 mg for IP, 20% of IP loaded) from HeLa cells. Antibodies: Affinity purified rabbit anti-TAF2 antibody used for IP at 10 mcg/mg lysate. TAF2 was also immunoprecipitated b... See more at [www.novusbio.com/NBP1-21372](http://www.novusbio.com/NBP1-21372)

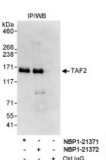

Updated 3/9/2013 2.0

**Limitations:** This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee).**

**Novus USA**

8100 Southpark Way, A-8  
 Littleton, CO 80120  
 p) 888-506-6887  
 p) 303-730-1950  
 f) 303-730-1966  
[novus@novusbio.com](mailto:novus@novusbio.com)

**Novus Canada**

461 North Service Road West,  
 Unit B37  
 Oakville, ON L6M 2V5  
 p) 855-668-8722  
 f) 905-827-6402  
[canada@novusbio.com](mailto:canada@novusbio.com)

**Novus Europe**

12 Cambridge Science Park  
 Cambridge, UK CB4 0FQ  
 UK: p) +44-1223-426-001, f) +44-871-971-1635  
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79  
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340  
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31  
[europe@novusbio.com](mailto:europe@novusbio.com)

