

NBP1-21405**TAB1 Antibody****Unit Size**

0.1 ml

Concentration

0.2 mg/ml

Storage

Store at 4C. Do not freeze.

Applications

WB, IP

Species

Human, Mouse

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.09% Sodium Azide

Immunogen

The immunogen recognized by this antibody maps to a region between residue 435 and 485 of human TAK1-binding protein 1 using the numbering given in entry NP_006107.1 (GeneID 10454).

Recommended Dilutions

Immunoprecipitation 2 -5 ug/mg lysate, Western Blot 1:2000 -1:10000

Buffer


Tris-buffered Saline and 0.1% BSA

Uses

This antibody is useful for western blot and immunoprecipitation.

Images

Western Blot: TAB1 Antibody [NBP1-21405] - Whole cell lysate from HeLa, 293T, and mouse NIH3T3 cells. NBP1-21405 used for WB at 0.04 mcg/ml (A) and 1 mcg/ml (B) and used for IP at 3 mcg/mg lysate.

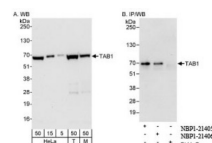

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

