

NB200-337**MED23 Antibody****Unit Size**

0.1 ml

Concentration

1.0 mg/ml

Storage

Store at 4°C. Do not freeze.

Applications

IP, WB (-)

Species

Human

Host

Rabbit

Clonality

Polyclonal

Purity

Immunogen affinity purified

Preservative

0.09% Sodium Azide

Specificity/Sensitivity

NB200-337 is specific for human CRSP 3 protein.

Immunogen

A synthetic peptide that maps to a region between residues 1 and 50 of human Cofactor Required for Sp1 transcriptional activation, subunit 3 using the numbering given in Swiss-Prot entry Q9ULK4 (GeneID 9439).

Recommended Dilutions

Immunoprecipitation 1 -4 ug/mg lysate
 Western Blot (Negative)

Buffer


Tris-citrate/phosphate buffer, pH 7 to 8

Uses

Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use. Western Blot: Not Recommended, use NB 200-339 for WB. *In some cases, the antibody may be diluted further than indicated.

Images

Immunoprecipitation: MED23 Antibody [NB200-337]
 - Nuclear extract (5 mg) from HeLa cells. Antibody used at 4 ug/5 mg extract.

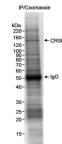

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

