

NB100-56393**TMED4 Antibody****Unit Size**

0.1 mg

Concentration

Please see the vial label for concentration.

Storage

Store at 4°C short term. Aliquot and store at -20°C long term. Avoid freeze-thaw cycles.

Applications

WB

Species

Human, Mouse, Rat

Host

Rabbit

Clonality

Polyclonal

Purity

Protein G purified

Preservative

0.05% Sodium Azide

Specificity/Sensitivity

An approx. 26 kDa protein is recognized in HL60 cell lysate.

Immunogen

A synthetic peptide containing amino acids [BLAST]CNYRTQMWDKQKEV[/BLAST] (aa55-67) of human HNLF. Genbank accession no. EAL23751.

Recommended Dilutions

Western Blot 1-3 µg/ml

Buffer


PBS containing 0.2% gelatin

Uses

Western Blot

Images

Western Blot: HNLF Antibody [NB100-56393] - Western blot analysis of HNLF in HL60 cell lysates using IMG-5112A at 2 µg/ml. Lane 1: Unabsorbed NB100-56393 Lane 2: Blocked with immunizing peptide

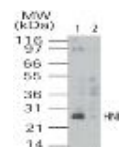

Updated 3/9/2013 2.0

Limitations: This product is for research use only and is not approved for use in humans or in clinical diagnosis. Products are guaranteed for six months from the date of receipt except for peptides and proteins which are guaranteed for 3 months. **For more information on our guarantee, please visit www.novusbio.com/guarantee.**

Novus USA

8100 Southpark Way, A-8
 Littleton, CO 80120
 p) 888-506-6887
 p) 303-730-1950
 f) 303-730-1966
novus@novusbio.com

Novus Canada

461 North Service Road West,
 Unit B37
 Oakville, ON L6M 2V5
 p) 855-668-8722
 f) 905-827-6402
canada@novusbio.com

Novus Europe

12 Cambridge Science Park
 Cambridge, UK CB4 0FQ
 UK: p) +44-1223-426-001, f) +44-871-971-1635
 DE: p) +49-6922-22340-60, f) +49-0800-58926-79
 IT: p) +39-02-4032-6786, f) +39-02-4032-6340
 FR: p) +33-1-76-77-45-30, f) +33-1-76-77-45-31
europe@novusbio.com

