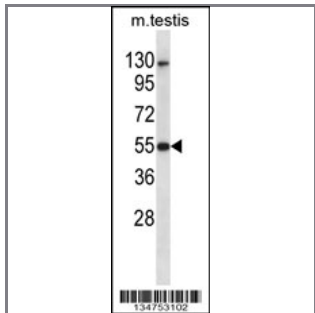


Mouse Dapk3 Antibody (Center)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)

Catalog #	Applications:	Reactivity:	Accessions:
AP14711c	WB, E	M	NP_031854.1 , NP_001177402.1 , O54784 , NP_001177403.1
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB34753
Application Data:	Calculated MW: 51422 Da		
	<p>Mouse Dapk3 Antibody (Center) (Cat. #AP14711c) western blot analysis in mouse testis tissue lysates (35ug/lane). This demonstrates the Dapk3 antibody detected the Dapk3 protein (arrow).</p>		
Gene ID:	Gene Symbol:		
13144	Dapk3		
Other Names:	Death-associated protein kinase 3;Dapk3;Zipk		
Target/Specificity:	<p>This Mouse Dapk3 antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 254-282 amino acids from the Central region of Mouse Dapk3.</p>		
Application Notes:	<p>The suggested dilution is: ELISA 1:1,000 Western blotting 1:100~500</p>		
Format:	<p>Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.</p>		

Storage:

Mouse Dapk3 Antibody (Center) can be refrigerated at 2-8°C for up to 6 months. For long term storage, place the at -20°C.

Precautions:

Mouse Dapk3 Antibody (Center) is for research use only and not for use in diagnostic or therapeutic procedures.
